




Janak **Product Catalog**

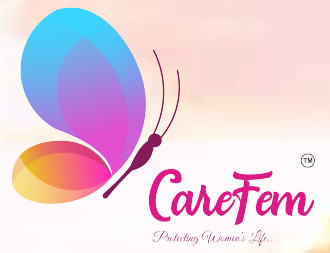


Flex Endosampler

Utilized for sampling of endometrial tissue, the unique design of the Flex Endosampler and allows for superior accuracy and a more reliable sample than traditional endometrial sample devices. The Endosampler is available with two holes for easier and more comfortable suction.




- Includes 3mm curette & 10cc syringe with lock-spring mechanism
- Highly-visible, graduated markings on curette
- Semi-rigid curette allows for easier insertion and injection of saline or contrast
- Curved design for easier manipulation
- Back-flow prevented by syringe lock-spring mechanism
- Sharper slot at curette tip for better sampling
- “Closed” system prevents sample contamination

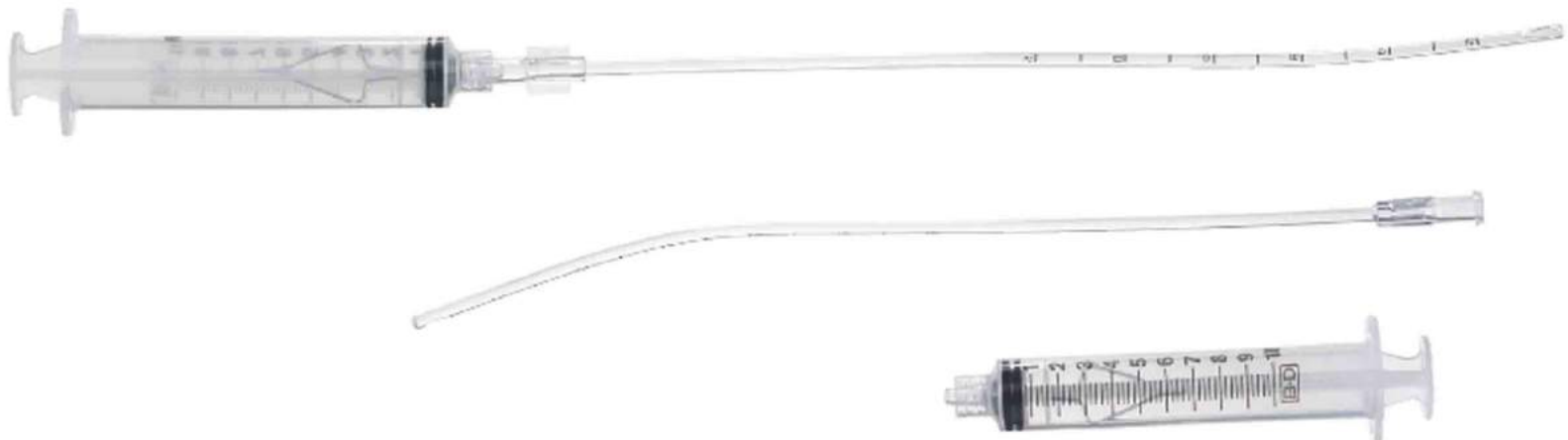

FlexPap Smear Devices


The Flex Healthcare pap smear device simultaneously samples both the endocervix and ectocervix eliminating the need for both a Cytobrush and a Spatula.

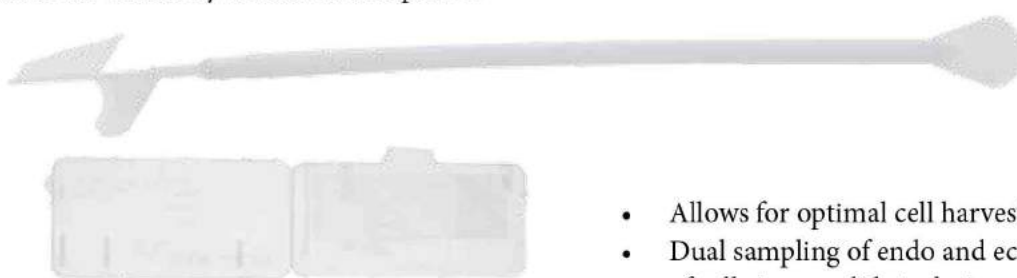

- Allows for optimal cell harvesting and cell transferring
- Dual sampling of endo and ectocervix and separation of cells in one-slide technique
- Minimizes false negatives
- Minimizes red blood cells and blood artifacts

Flex Cytobrush and Pap Cell

Standard devices for one of the most frequently performed gynecology procedures.

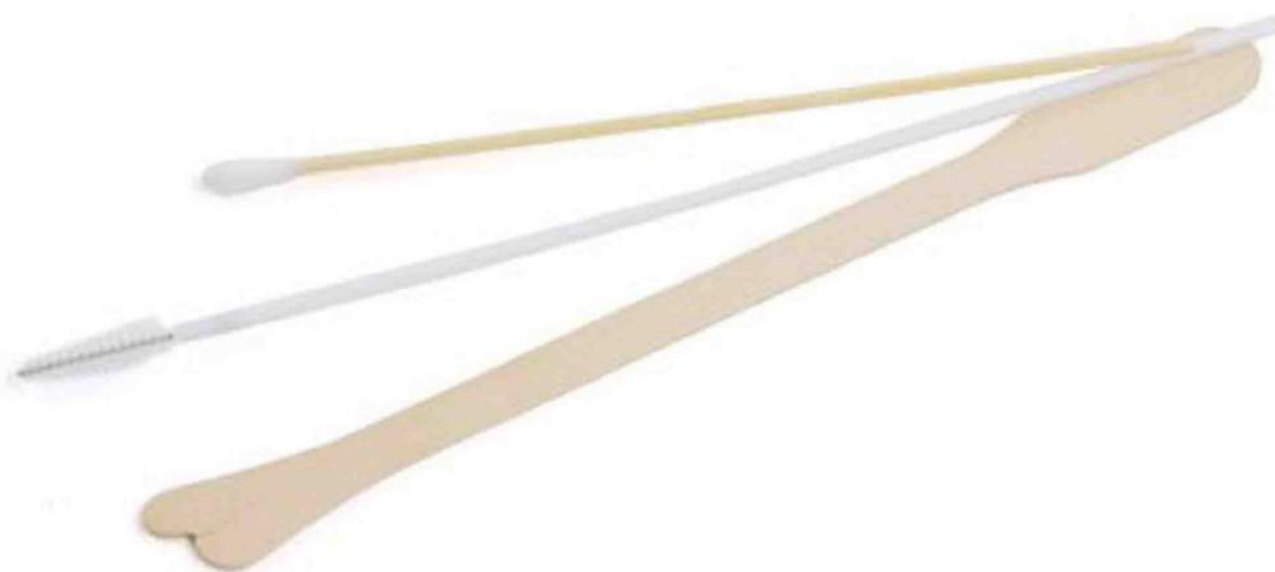

Flex IUI Catheters

Flex Insemination Catheter is a sterile, flexible catheter designed for use in intrauterine insemination. Tip of catheter is rounded and beveled for patient comfort.



Regular




Mini



Curved

Flex HSG Catheters

Flex HSG Catheter is used to inject contrast media or saline into the uterus and fallopian tubes for hysterosalpingography and sonohysterography procedures. It has a latex-free balloon catheter with a rigid insertion sheath.

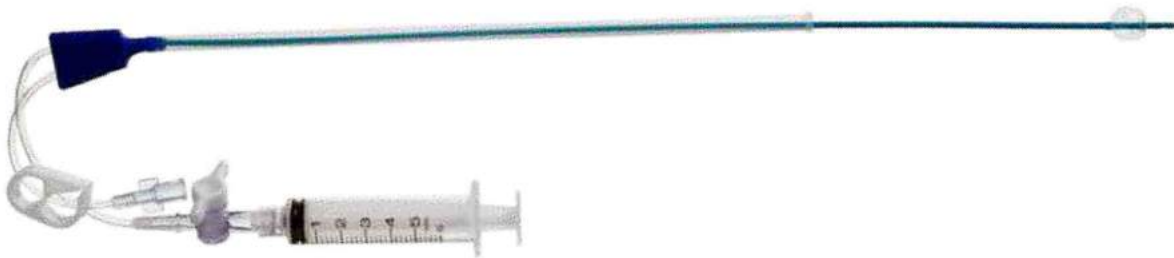

Flex LED Speculum

Flex LED Speculum has a cordless LED light to improve visualization of the exam area.



Large




Medium



Small

Flex Disposable Electrodes

Made of high-quality Tungsten wire, for superior shape and integrity, Flex offers a full range of shapes and sizes to suit most electrosurgical procedures.

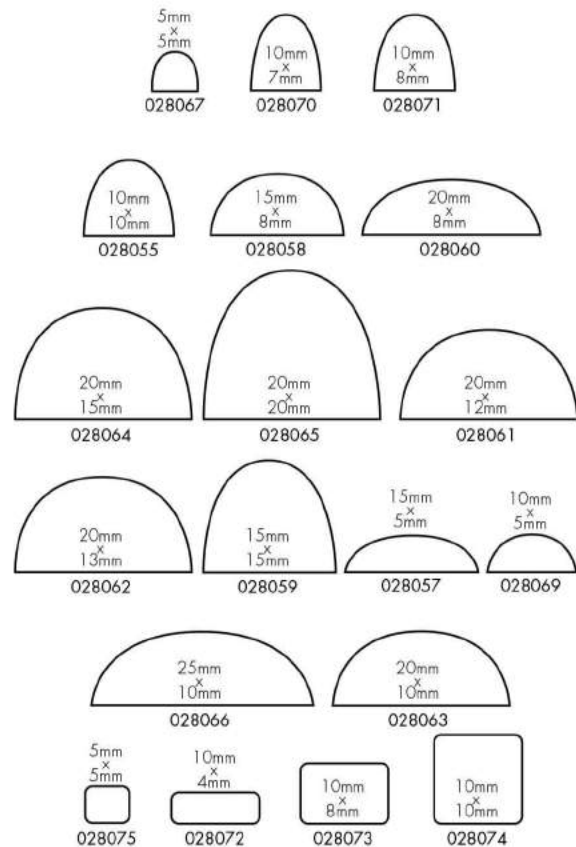

Flex Disposable Hand Switch Pencil



Vaginal Speculum

Disposable Speculum is available in both Graves and Pederson styles and easily locks in a variety of positions to reduce examination time.



Large



Medium



Small

Vacuum Assisted Delivery Cups

Vacuum Delivery Cups offer optimal adhesion to ensure proper vacuum levels



Release Valve

Umbilical Cord Clamp

Umbilical Cord Clamp has grooved clamping and lock to secure the umbilical cord prior to cutting.




Ipas

Manual Vacuum Aspiration

AN EFFECTIVE, SAFE, PATIENT-CENTERED SOLUTION THAT
USES SUCTION FOR UTERINE EVACUATION

IT IS A COST-EFFECTIVE METHOD THAT DOES NOT
REQUIRE ELECTRICITY, FEASIBLE IN OUTPATIENT
SETTINGS AND CAN BE DONE UP TO 12 WEEKS


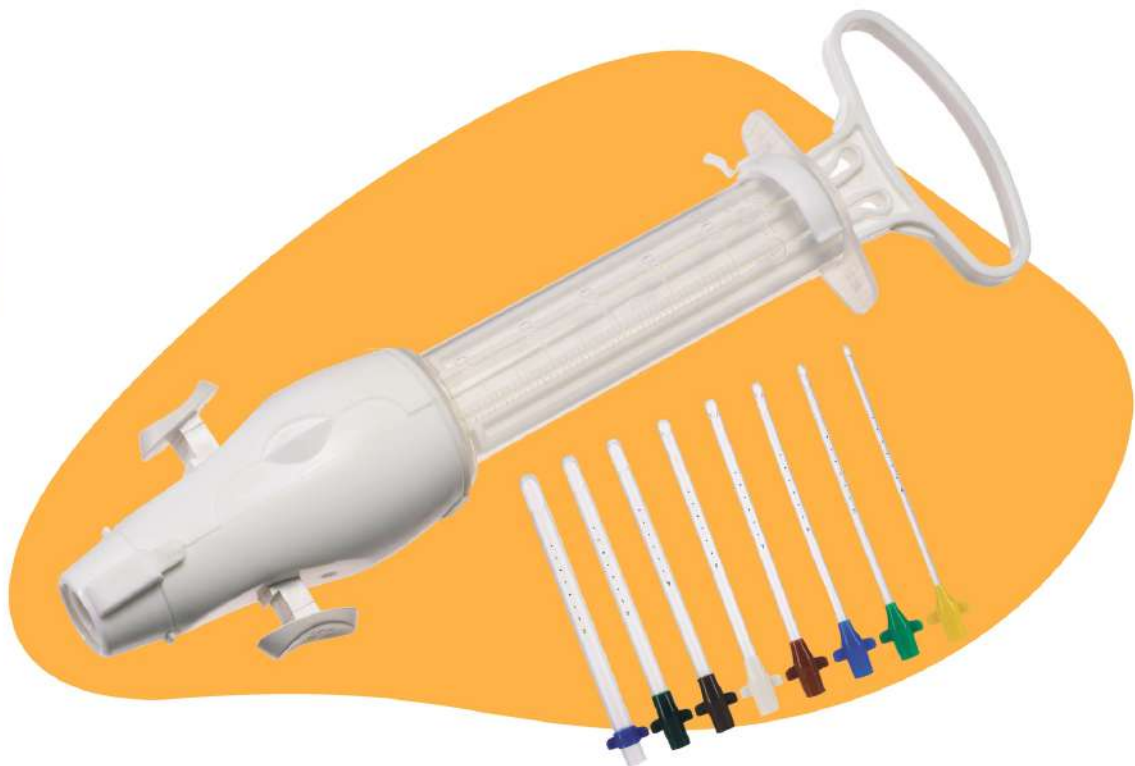


WomanCare GLOBAL | A COMPANY OF
dkt INTERNATIONAL

Flexible/Karman Cannula-Curette

- 022004K 4mm (10/Pack)
- 022005K 5mm (10/Pack)
- 022006K 6mm (10/Pack)
- 022007K 7mm (10/Pack)
- 022008K 8mm (10/Pack)
- 022009K 9mm (10/Pack)
- 022010K 10mm (10/Pack)
- 022011K 11mm (10/Pack)



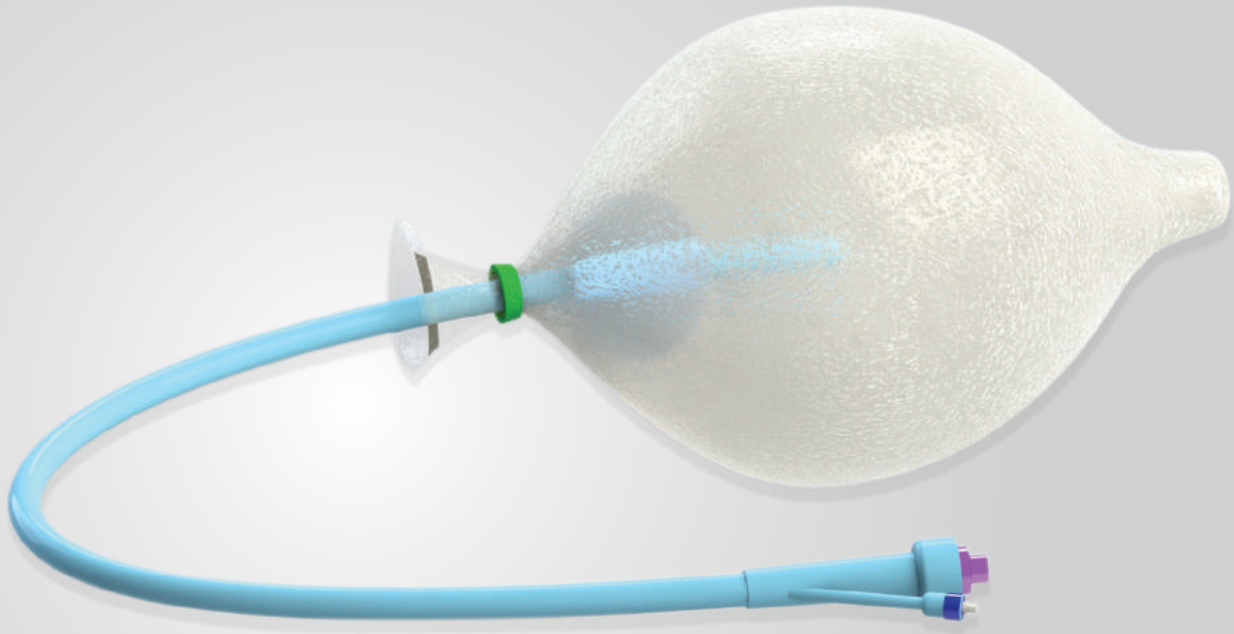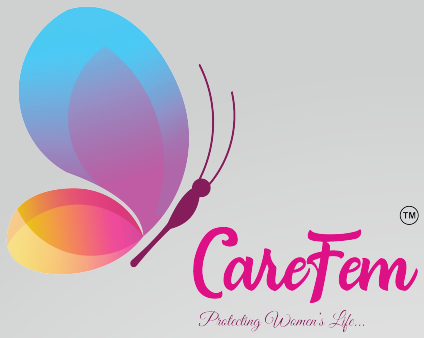

Dr. Burke's ESM-UBT Kit
UTERINE BALLOON TAMPONADE

Pregna is proud to team up with Vayu Global Health Foundation to introduce Dr. Burke's ESM UBT Kit

Features:



At a time when Every Second Matters, Dr. Burke's ESM-UBT Kit eliminates the time required in procuring the required components.



It reduces the expense of assembling, procuring and reprocessing the required components from various parties.



Single use device reduces the risk of cross-contamination and reprocessing.



Effective in rapid treatment of PPH, within 5 to 15 minutes after insertion and inflation.



Each kit is assembled with all the components required to adequately handle PPH, making it suitable even for low-resource settings



With all the components packed in a easy to store kit, it increases the efficiency of medical staff in the event of an emergency




The compactness of the kit allows it to be easily carried from one location to another




Easy to assemble with minimal training required.

Dr. Burke's ESM-UBT Kit is a cost effective and widely available solution which helps maximize its reach to every corner of the country.


DR. BURKE'S ESM-UBT KIT PLACEMENT

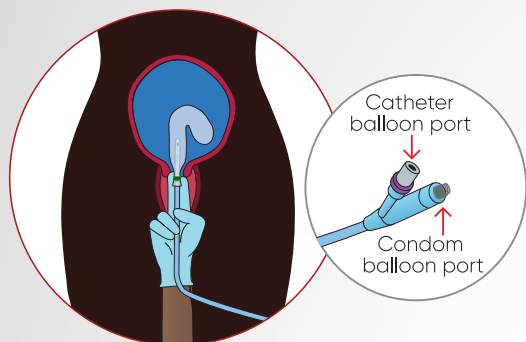


STEP 1: Locate the cervix and insert the device.

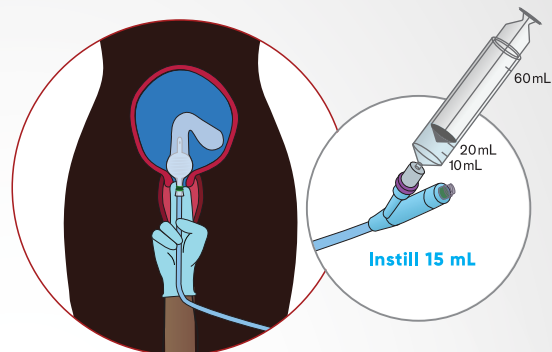

STEP 2: Inflate the balloon of the urinary catheter with saline.

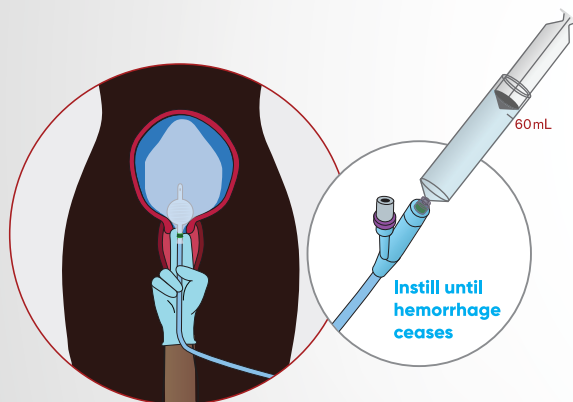

STEP 3: Inflate the condom balloon with saline.
The maximum balloon volume is 500 ml.

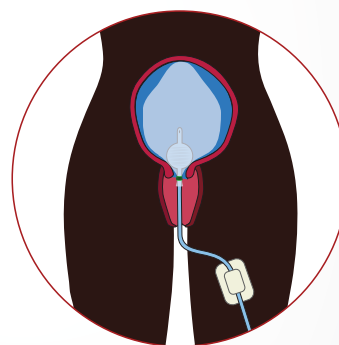

STEP 4: Secure with the catheter holder.

Download Vayu UBT Mobile App for more information



Enroll for free: ESM-UBT device training course for medical practitioners



eloira

Hormonal Intrauterine System
Levonorgestrel 52 mg



Eloira is India's first indigenous progestin-only reversible intrauterine contraceptive system having a 99.8% success rate and is effective upto 5 years.

Eloira consists of a T-shaped plastic frame with a reservoir containing the drug Levonorgestrel (LNG) which sits in the uterine cavity.


The reservoir is provided with a rate controlling membrane that releases a minimal amount of LNG (approx. 20 micrograms per day) into the uterus.

Contains 52 mg Levonorgestrel with release of approx. 20 µg per day



How does it work?

Eloira works in 3 ways

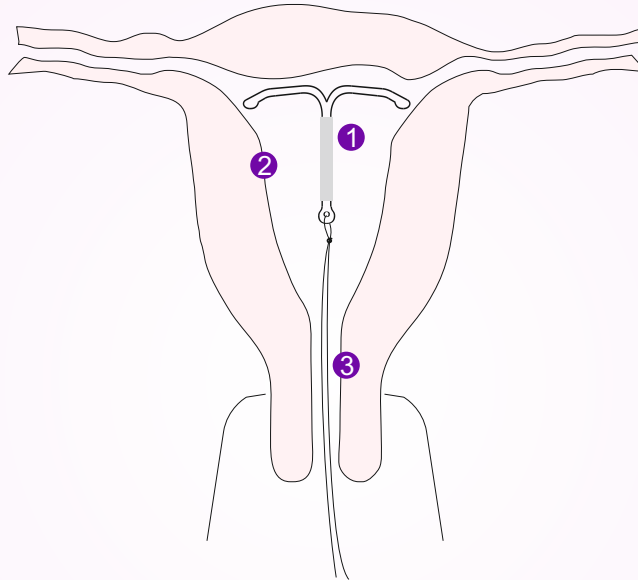

- 1 A weak local foreign body effect makes it an unfavorable environment in the uterus to accept an egg
- 2 Preventing build-up of the endometrium lining
- 3 Thickens the mucus lining in the cervix thereby inhibiting mobility of the sperms.

Eloira comes with several non-contraceptive benefits too.

Highly effective upto 99.8%

Reversible

Effective for upto 5 years

No long term side effects

Lighter periods and less blood loss

Reduced menstrual cramps

Lowers the risk of anemia

Affordable price and readily available



CryoPop [®]


An avant garde cryotherapy device

- A Cryotherapy device for treatment of precancerous cervical lesions
- No cables required during operation
- Solid CO₂ based system consumes upto 75% lesser gas
- Solid state conduction makes it highly cooling efficient



Dr. Burke's ESM-UBT Kit

- Uterine Balloon Tamponade
- Effective in rapid treatment of Postpartum Bleeding
- Single use device reduces the risk of cross contamination
- Easy to assemble with minimal training required
- Majority of deaths are preventable

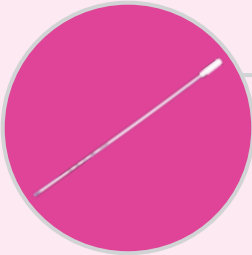

Accurette [®]

Endometrial Biopsy Curette

- Pain free collection of endometrial tissue for pathology.
- Safe and simple procedure.
- Sufficient sample quantity, eliminate need of repetition.



esā [®]

Disposable Uterine Sound

- Smooth finish with special rounded tip to avoid perforation of fundus.
- Wide handle with serrations to give better grip
- Single use device to avoid cross contamination.
- Anatomical design with high flexibility to suit the uterine cavity

Pregna's Range Of IUDs

All of our products are saliently designed with:

- High quality soft nylon suture to provide utmost comfort to the patient.
- Medical grade LDPE frames blended with Barium Sulphate for visibility under X-ray.
- Depth gauge marking on insertion tube to provide accuracy and convenience while insertion.
- Reversible contraception with 99% effectivity



eloira® 52 mg Levonorgestrel (LNG IUS)

- Effective period: 5 years
- Constant release of 20mg of LNG per day
- Flexible arms for perfect fundal placement
- Used in treatment of HMB (AUB & DUB) and fibroids less than (<3cm).
- Comes with Esa - disposable uterine sound



etherena® T Cu 380A

- Effective period: 10 years
- Specially designed Intrauterine Enabler (IUE) for aseptic loading
- Single handed operation
- Precise uterine depth adjuster for added accuracy
- Comes with Esa - disposable uterine sound.



PREGNA+® T Cu 380A

- Effective period: 10 years
- Simple and safe insertion & removal
- Preloaded solid rod for easy insertion.



Inara[®] Cu 375 / Cu 375 Sleek

Effective period: 5 years
Specially designed inverted 'U' shape frame for added comfort while inserting.
Preloaded and ready for insertion
Inara Cu 375 Sleek has a specially designed smaller size 'U' shape frame
Inara Cu 375 Sleek is recommended for nulliparous women.
Both models come with Esa - disposable uterine sound in pack.



Inara[®] Cu 250 / Cu 250 Sleek

Effective period: 3 years
Specially designed inverted 'U' shape frame for added comfort while inserting.
Preloaded and ready for insertion
Low copper content
Inara Cu 250 Sleek has a specially designed smaller size 'U' shape frame
Inara Cu 250 Sleek is recommended for nulliparous women.
Both models come with Esa - disposable uterine sound in pack.


Post - Partum Range of IUDs


Specially designed preloaded long inserter for easy fundal placement with no touch technique.
Long nylon suture which is visible immediately after the insertion.
To be used within 48 hours of normal delivery.




 Cu 375

Effective period: 5 years

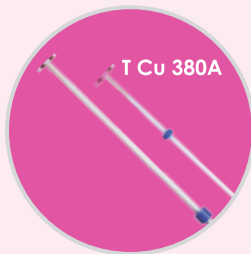

 T Cu 380A

Effective period: 10 years

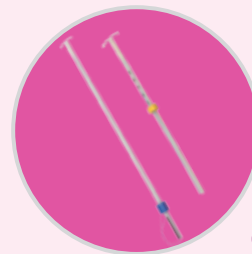

eloira 

Effective Period: 5 years

MOTHER NEEDS

Advanced Manual Breast Pump Premium Quality
Model No.: BB7069



HAPPY BREAST FEEDING

Small and easy to handle



Silicon pad for more comfort



Adjustable suction strength



66

MOTHER NEEDS



Breastfeeding Smile Cycle

With BuddsBuddy Breast pumps you can easily have a simple and gentle pumping. This product helps to easily pump and store the breast milk for your baby. Breast Pumps

The breast pump is especially made for the Mothers who are busy with work and are unable to be present in home. They can benefit their babies with their breast milk by using these breast pumps.

3 Key steps to effective breast milk expression

1. Fit

BuddsBuddy Breast pumps comes with air-tight silicone ring which can be easily adjustable according to the breast shape.

2. Gentle stimulation

Gentle stimulation reduces the burden on the Mother's breasts. The pump processes with low suction mode.

3. Expression

While pumping the milk is collected in the container safely. You can control the pump and easily adjust while using it

How to use a Breast Pump?

When Mother's special hormones are stimulated, breast milk is produced. When the baby is breastfed the hormones are released. When using the breast pump, the Mother needs to be relaxed, thus the pump will promote lactation. Stop pump in case you are feeling pain.

COMFORTABLE



Relaxed and comfortable position for an effective milk expression.

APPROVED EFFECTIVENESS



Adjustable expression, ergonomic handle.

SIMPLE

CONVENIENT

CLEVER

67

MOTHER NEEDS

Premium Silicone Compact Breast Pump 1pc
Model No.: BB7073
100% SILICONE MATERIAL




Silicone Nipple Puller With Case 1pc
Model No.: BB7067
100% SILICONE MATERIAL

Use BuddsBuddy Nipple Puller to gently draw your nipples out to breastfeed the baby. The product is absolutely easy to clean and use and don't have any side effects as well.

How to use

Wash the nipple puller before using it. Place the shield around the nipple. Press the bulb and push the shield down. Repeat until the nipple is fully extended

The effect of after use

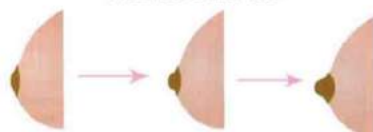

Safe, Convenient

Hassle Free Process

Easy To Use

68

MOTHER NEEDS

Premium Comfort Nursing Breast Shell 2pcs
Model No.: BB7068

Happy breast feeding to you



Soft Silicone Protective Nipple Shield With Case 2pcs
Model No.: BB70074

Enjoy the Happiness of breast feeding

BuddsBuddy Nipple shield is made from soft silicone. Imitates similar like Mother's natural nipples.

How to use

Place the shield over the nipples on breasts. Hold the shield with fingers while feeding. Breastfeed frequently for short periods. Before touching the nipples, wash the hands thoroughly. Warm the breasts with hot towel to stimulate breast milk. Massage the nipples to open mammary glands before using a nipple shield.

1. Clean the nipple first



2. Drop a little breast milk into the nipple



3. Press it from the left and right sides and press it from above



69



MOTHER NEEDS

Honeycomb Pad

This special lining allows the excess milk to spread evenly. The leak-proof cover protects your clothing, thus keeps your breasts dry for a longer time.

Quick Absorbency, Breathable layer, Leak-proof coverage

Breast pads are quick absorbent and contains breathable layer. Ultra slim in size in fits easily under the clothing. Absolutely leak-proof, can keep you dry for hours.

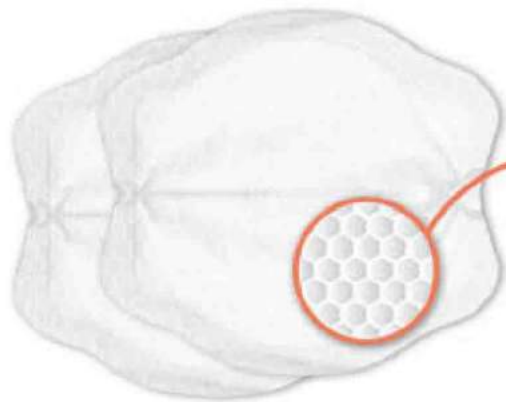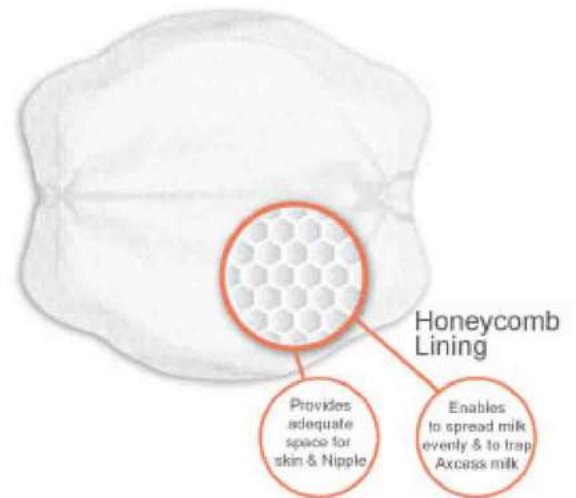

HONEYCOMB DESIGN NEW BREATHABLE LINING

Why BuddsBuddy?

- * FABRIC WITH ANTIBACTERIAL
- * SUPER ABSORBENT
- * MAXIMUM BREATHABILITY
- * THIN AND INVISIBLE



COTTON SOFT



QUILTED HONEYCOMB LINING



OVERNIGHT PROTECTION

HOW TO USE



1
Removal pad from individual packaging.




2
Peel off the protective backing from the nonslip tape.



3
Face the nonslip tape side against bra cup.



4
Position and press pad against the bra cup.

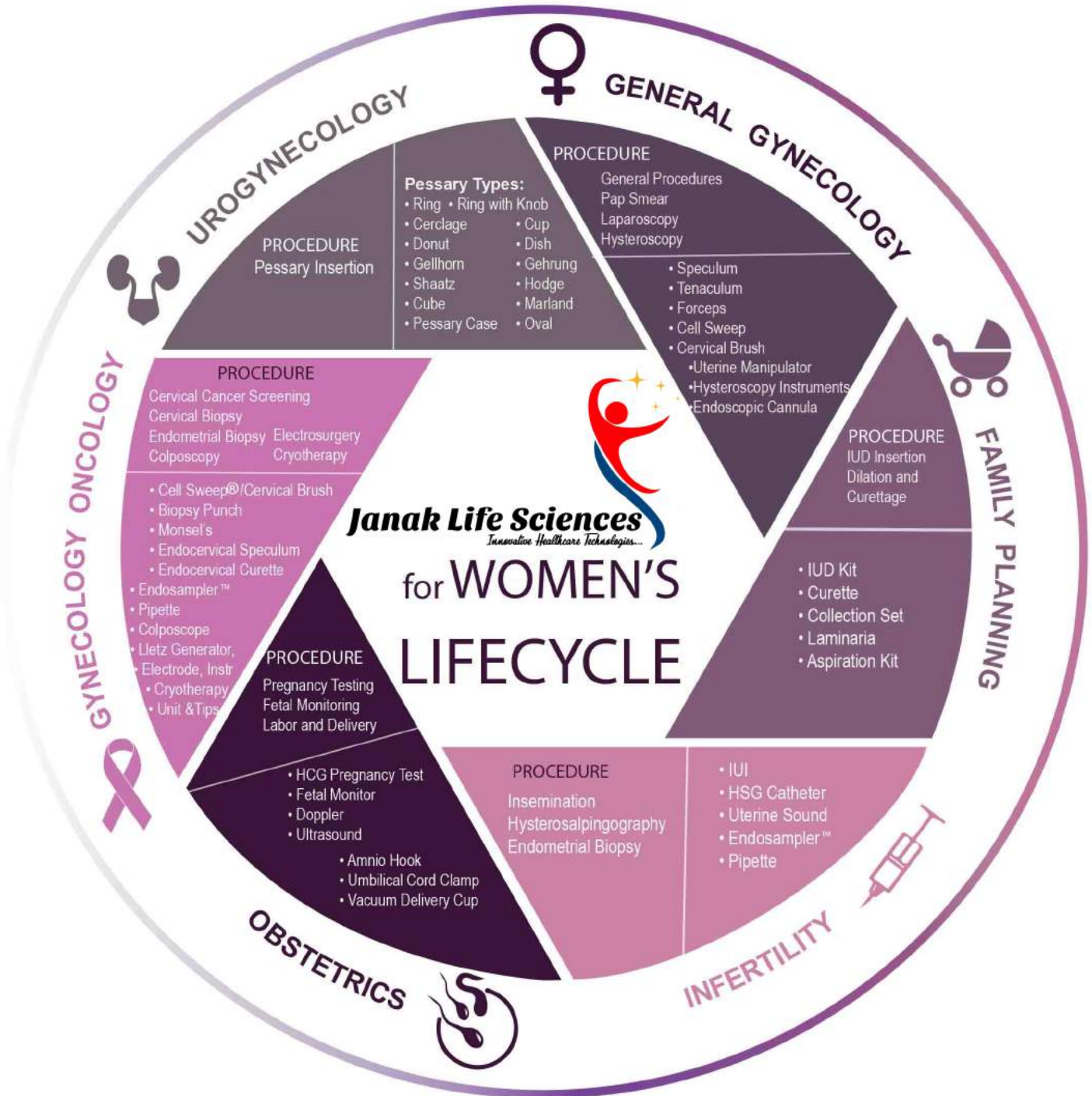

UROGYNECOLOGY

PROCEDURE
Pessary Insertion

- Pessary Types:**
- Ring
 - Ring with Knob
 - Cerclage
 - Donut
 - Gellhorn
 - Shaatz
 - Cube
 - Pessary Case
 - Cup
 - Dish
 - Gehrung
 - Hodge
 - Marland
 - Oval

GENERAL GYNECOLOGY

PROCEDURE
General Procedures
Pap Smear
Laparoscopy
Hysteroscopy

- Speculum
- Tenaculum
- Forceps
- Cell Sweep
- Cervical Brush
- Uterine Manipulator
- Hysteroscopy Instruments
- Endoscopic Cannula

FAMILY PLANNING

PROCEDURE
IUD Insertion
Dilation and Curettage

- IUD Kit
- Curette
- Collection Set
- Laminaria
- Aspiration Kit

GYNECOLOGY ONCOLOGY

PROCEDURE
Cervical Cancer Screening
Cervical Biopsy
Endometrial Biopsy
Colposcopy

- Cell Sweep@/Cervical Brush
- Biopsy Punch
- Monsel's
- Endocervical Speculum
- Endocervical Curette
- Endosampler™
- Pipette
- Colposcope
- Lletz Generator
- Electrode, Instr
- Cryotherapy
- Unit & Tips

PROCEDURE
Pregnancy Testing
Fetal Monitoring
Labor and Delivery

- HCG Pregnancy Test
- Fetal Monitor
- Doppler
- Ultrasound

- Amnio Hook
- Umbilical Cord Clamp
- Vacuum Delivery Cup

OBSTETRICS

PROCEDURE
Insemination
Hysterosalpingography
Endometrial Biopsy

- IUI
- HSG Catheter
- Uterine Sound
- Endosampler™
- Pipette

INFERTILITY